

ALBAhance™ LISS ADDITIVE REAGENT

REF Z333

- No US Standard of Potency
- Discard if turbid
- Preservative: 0.09% sodium azide

CAUTION: THIS PRODUCT HAS COMPONENTS (DROPPER BULBS) CONTAINING DRY NATURAL RUBBER


INTERPRETATION OF LABELING SYMBOLS



Batch code



Use by (YYYY-MM-DD)

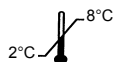

Storage temperature limitation (2°C-8°C)



In vitro diagnostic medical device



Consult instructions for use



Product code



Manufacturer

INTENDED USE

ALBAhance™ LISS Additive Reagent is intended for use as a potentiator in antibody detection, antibody identification and compatibility test procedures.

SUMMARY AND EXPLANATION

A substantial reduction in the incubation time for antigen/antibody mixtures can be achieved when the red blood

cells and serum are suspended in a LISS medium. It is also recognized that most antibodies will show an increase in the test sensitivity when LISS is incorporated into the test medium.

PRINCIPLE OF THE TEST

ALBAhance™ LISS Additive Reagent is added directly to antibody detection, antibody identification or cross-match reagents to reduce the ionic strength of the testing environment. There is enhancement of antigen-antibody interactions during incubation. Because antibody uptake is enhanced, incubation periods of low ionic test systems are generally shorter than those of routine saline/albumin tests.

REAGENT DESCRIPTION

ALBAhance™ LISS Additive Reagent is a low ionic reagent containing 0.09% sodium azide w/v as a preservative. ALBAhance™ LISS Additive Reagent is to be used as supplied following the directions detailed in this insert. The volume delivered by the reagent dropper bottle is approximately 40µL. Bearing this in mind, care should be taken to ensure that appropriate serum:cell ratios are maintained in all test systems.

PRECAUTIONS

Store at 2°C - 8°C.
Do not use if turbid.
Do not dilute.

Do not use beyond expiry date.
This reagent contains 0.09% (w/v) sodium azide. Sodium azide may be toxic if ingested and may react with lead and copper plumbing to form explosive compounds. If discarded into sink, flush with a large volume of water to prevent azide build-up.

This product has components (dropper bulbs) containing dry natural rubber.

This reagent is for *in vitro* diagnostic use only.

SPECIMEN COLLECTION AND PREPARATION

Specimens should be collected by aseptic technique with or without an anticoagulant. The specimen should be tested as soon as possible after collection. If testing is delayed, the specimen should be stored at 2°C - 8°C. Blood specimens exhibiting gross hemolysis or contamination should not be used. Do not use collection tubes that contain serum or plasma/cell separation media. Clotted samples or those collected in EDTA should be tested within fourteen days from collection. Donor blood stored in citrate anticoagulant may be tested until the expiry date of the donation.

MATERIALS

Materials provided

- ALBAhance™ LISS Additive Reagent

Materials required but not provided

- Patient/donor serum/plasma
- Isotonic saline
- Reagent red blood cells for antibody detection or identification
- Antiglobulin Reagent
- IgG sensitized red blood cells
- 10 x 75mm or 12 x 75mm glass test tubes
- Pipettes
- Centrifuge
- Heating block / waterbath @ 37°C
- Timer
- Optical aid

TEST PROCEDURE

General Information

This reagent has been standardized for use by the technique described below and therefore its suitability for use in other techniques cannot be guaranteed. When a test is required to be incubated for a specific period of time, a timer should be used.

RECOMMENDED TECHNIQUE

37°C Indirect Antiglobulin Test

- Add 2 drops of patient's serum/plasma to a test tube.
 - Add 1 drop of red blood cells suspended to 2-4% in isotonic saline.
 - Add 2 drops of ALBAhance™ LISS Additive Reagent
 - Mix the test well and incubate for 15-45 minutes at 36-38°C.
 - Wash the test at least 3 times with a large excess isotonic saline. e.g. 4ml of saline per 12 (or 10) x 75mm glass tube)
- NOTE:** (i) allow adequate spin time to sediment the red blood cells.
(ii) make sure that most of the residual saline is removed at the end of each wash.
- Add Antiglobulin Reagent to each test tube in the amount specified in the manufacturer's product insert.
 - Mix the contents of the test tube well and centrifuge. Suggested centrifugation: 900-1000g for 10 seconds or a time and speed appropriate for the centrifuge used that produces the strongest reaction of antibody with antigen-positive cells, yet allows easy resuspension of antigen-negative cells.
 - Gently shake the test tube to dislodge the cell button from the bottom and observe macroscopically for agglutination. Negative reactions may be examined with an optical aid.
 - Record results.
 - Add IgG sensitized antiglobulin control cells to confirm the validity of negative test results.

STABILITY OF REACTION

Test results should be read and interpreted immediately after centrifugation. Delays may cause dissociation of antigen-

antibody complexes resulting in weak positive or false negative reactions.

INTERPRETATION OF RESULTS

Agglutination = positive test result
No agglutination = negative test result

QUALITY CONTROL

Quality control of reagents is essential and should be performed on each day of use and in accordance with local, state and federal regulations.

PERFORMANCE LIMITATIONS

Direct antiglobulin test positive samples will react by the indirect antiglobulin test irrespective of their antigen status.

Driblocks and waterbaths promote better heat transfer and are recommended for 37°C tests, particularly where the incubation period is 30 minutes or less.

Gently resuspend tube tests before reading. Excessive agitation may disrupt weak agglutination and produce false negative results.

Excessive centrifugation can lead to difficulty in resuspending the cell button, while inadequate centrifugation may result in agglutinates that are easily dispersed.

The expression of certain red blood cell antigens may diminish in strength during storage, particularly in EDTA and clotted samples. Better results will be obtained with fresh samples.

False positive or false negative results can occur due to contamination of test materials, improper reaction temperature, improper storage of materials, omission of test reagents and certain disease states.

SPECIFIC PERFORMANCE CHARACTERISTICS

Prior to release, each lot of ALBAhance™ LISS Additive Reagent is tested by the method detailed in the package insert against a panel of antigen-positive and antigen-negative red blood cells to ensure suitable reactivity.

For additional information or technical support, contact Product Technical Support at 1-888-228-1990.

BIBLIOGRAPHY

1. Low B, Messeter L. Antiglobulin test in low ionic strength salt solution for rapid antibody screening and crossmatching. Vox Sang 1974; 26:53.
2. Moore HC, Mollison PL. Use of a low ionic strength medium in manual tests for antibody detection. Transfusion 1976; 16:291.

DATE OF ISSUE

July 2010

US Distributor:

Quotient Biodiagnostics Inc.
41 University Drive
Newtown
PA 18940
USA

Customer Service Tel: 1-888-284-1901
Product Technical Support Tel: 1-888-228-1990
Customer Service Fax: 1-888-694-5208
E-Mail: customer.service@quotientbd.com
Web: www.quotientbd.com

Manufacturer:

Alba Bioscience Limited
Ellen's Glen Road
Edinburgh
Scotland, UK
EH17 7QT

Tel: +44 (0) 131 536 5907
Fax: +44 (0) 131 536 5897
E-Mail: customer.services@albabioscience.co.uk

© Alba Bioscience Limited 2010

Z333QPI/02