



# QUOTIENT

## ALBAcyte®

### B Cells

#### REAGENT RED CELLS FOR REVERSE GROUPING

REF Z411



1434


#### INTERPRETATION OF LABEL SYMBOLS



Batch code



Use by (YYYY-MM-DD)

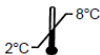

Storage temperature limitation (2-8 °C)



*In vitro* diagnostic medical device



Product code



Consult instructions for use

[www.quotientbd.com](http://www.quotientbd.com)



Manufacturer

#### INTENDED PURPOSE

These reagent red cells are for the ABO reverse grouping of patient or donor serum/plasma.

#### INTRODUCTION

ABO blood grouping is generally performed by testing red cells with anti-A and anti-B (many laboratories also test with anti-A,B). A<sub>1</sub>, A<sub>2</sub>, B and O cells should be used to control each batch of tests and single tests. Confirmation of the red cell group is normally provided by simultaneously performing a reverse or serum group ie testing the donor or recipients serum/plasma with reagent red cells of groups A<sub>1</sub> and B.

#### REAGENT DESCRIPTION

These reagent red cells are presented as a 2-3% suspension of pooled washed red cells in Modified Alsever's Solution. The Rh phenotype of the B red cells is ccdee. The preservative solution has been specially formulated to preserve red cell integrity and antigenicity and contains the following components - trisodium citrate, citric acid, dextrose, inosine,

neomycin sulphate (0.103g/l) and chloramphenicol (0.349g/l). These reagent red cells may be used directly from the vial or may be washed and resuspended before use to 2-3% in PBS or 1.5-2% in LISS. Reagent red cells treated in this way must be discarded within 24 hours of preparation. Transfer of these reagent red cells to another container is not recommended. Furthermore, when the user changes the reagent in any way, eg the preparation of LISS cell suspensions, the user is responsible for assuring the strength of red cell suspension, the quality of PBS or LISS used and the generation and storage of relevant documentation.

The volume delivered by the reagent dropper bottle is approximately 40µl; bearing this in mind, care should be taken to ensure that appropriate serum: cell ratios are maintained in all test systems.

**This reagent complies with the requirements of Directive 98/79/EC on *in vitro* Diagnostic Medical Devices and the recommendations contained in the Guidelines for Blood Transfusion Services in the United Kingdom.**

#### PRECAUTIONS FOR USE AND DISPOSAL

Source material from which this product is derived was found non reactive for HBsAg, Anti-HIV 1/2 and Anti-HCV.

**No known test method can offer assurances that products derived from human blood will not transmit infectious disease, therefore appropriate care should be taken in the use and disposal of this product.**

Chloramphenicol is classified as a carcinogen and neomycin sulphate is classified as an irritant.

This reagent is for *in vitro* professional use only.

#### STORAGE CONDITIONS

The reagent should be stored at 2°C - 8°C. Do not freeze. Do not use if obviously discoloured or haemolysed. Do not use beyond the notified expiry date.

#### SPECIMEN COLLECTION AND PREPARATION

Specimens should be collected by aseptic technique. The specimen should be tested as soon as possible after collection. If testing is delayed, the specimen should be stored at 2°C-8°C for a maximum of 48 hours. Blood specimens exhibiting gross haemolysis or contamination should not be used.

#### TEST PROCEDURES

No specific test procedures are recommended. Users are advised to carefully validate procedures and confirm reagent suitability before use.

#### PERFORMANCE LIMITATIONS

The presence of irregular antibodies in the serum/plasma of a patient/donor may cause unexpected agglutination of these reagent red cells.

For samples showing discrepant results, the patient's serum/plasma should be retested with their own red cells (autotest) and with group O red cells at room temperature. Some loss of antigenic expression may occur during the stated shelf life. Since this loss is partly determined by characteristics of individual blood donations or donors which

cannot be predicted or controlled, the recommended conditions of storage and use must be rigidly applied.

False positive or false negative results can occur due to contamination of test materials, improper reaction temperature, improper storage of materials, omission of test reagents and certain disease states.

#### SPECIFIC PERFORMANCE CHARACTERISTICS

The reagent red cells have been shown to have a negative direct antiglobulin test, indicating that no human IgG or C3 complement components are detectable on the cell surface.

#### DATE OF ISSUE

2023-08

For further information or advice please contact your local distributor.



**Emergo Europe B.V.**  
Westervoordijk 60  
6827 AT, Arnhem  
The Netherlands



**Alba Bioscience Limited**  
James Hamilton Way  
Penicuik  
EH26 0BF  
UK

Tel No: +44 (0) 131 357 3333  
Fax No: +44 (0) 131 445 7125  
E-Mail: [customer.servicesEU@quotientbd.com](mailto:customer.servicesEU@quotientbd.com)

© Alba Bioscience 2023

Z411PI/08

#### INSTRUCTIONS FOR USE

Instructions for Use are available on Quotient website at [www.quotientbd.com](http://www.quotientbd.com) or can be requested from your Local Distributor, by providing the relevant Product Code stated on product labels and Instructions for Use supplied with the product.

#### ISTRUZIONI PER L'USO

Le istruzioni per l'uso sono disponibili sul sito Web di Quotient, [www.quotientbd.com](http://www.quotientbd.com), o possono essere richieste al proprio distributore locale fornendo il codice del prodotto indicato sulle etichette del prodotto e nelle istruzioni per l'uso fornite con il prodotto.

#### GEBRAUCHSANWEISUNG

Die Gebrauchsanweisung ist auch auf der Quotient Internetseite unter [www.eu.quotientbd.com/ifu](http://www.eu.quotientbd.com/ifu) erhältlich oder kann bei Ihrem zuständigen Vertriebspartner angefordert werden. Hierfür geben Sie bitte die jeweilige Artikel-Nummer an, die sich auf den Etiketten und Gebrauchsanweisungen befindet, die mit dem Produkt geliefert werden.

#### ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ

Οι οδηγίες χρήσης είναι διαθέσιμες στον ιστότοπο της Quotient στη διεύθυνση [www.quotientbd.com](http://www.quotientbd.com) ή μπορείτε να τις ζητήσετε από τον τοπικό διανομέα, παρέχοντας τον αντίστοιχο κωδικό προϊόντος που αναφέρεται στις ετικέτες του προϊόντος και στις οδηγίες χρήσης που συνοδεύουν το προϊόν.